


PowerSpin™ BX Centrifuge

Product Images



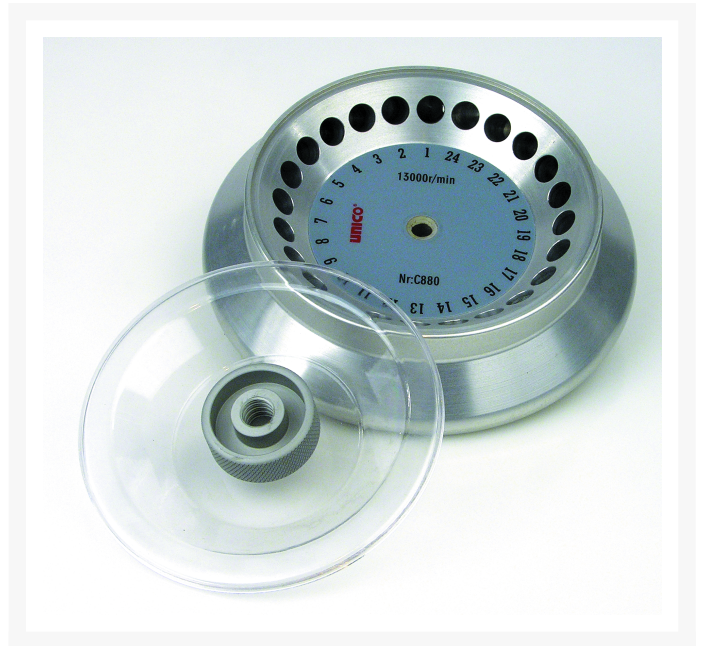

Short Description

1. POWERSPIN™ BX CENTRIFUGE
2. VARIABLE SPEED 1000-13000 RPM WITH 8 PLACE FIXED ANGLE, 24 PLACE MICROHEMATOCRIT AND 24 PLACE MICROCENTRIFUGE ROTORS
3. 99 MIN DIGITAL TIMER
4. POWER INPUT 110 V

Description

UNICO is proud to introduce the latest addition to the PowerSpin™ Family – the BX model high speed centrifuge with 3 rotor options:

8-place Blood Tube rotor, 24 place Microhematocrit and 24 place Micro-Tube rotor

The state of the art design incorporates advanced safety features like Zero RPM Lid Lock and Auto Unbalance Detection, Rotor Recognition and Shutdown.

The built in oversized LCD Display provides real time data to the user. With the push of a button you can speed up or slow down the rotor within 1,000 - 11,000 RPM range (Max 5000 RPM with the Blood Tube Rotor). Maintenance free brushless DC motor is powerful and super quiet. The cover of the Microhematocrit rotor incorporates a reader, allowing for instant reading and eliminating the need for additional specimen handling.

| | | |
|-----------------|-----------------|---|
| Accessories | Header | Accessories |
| | C880-08 | Fixed Angle Rotor, 8 place |
| | C880-02 | Microhematocrit Rotor, 24 place |
| | C880-03 | Microtube rotor, 24 place |
| | C880-10 | BX Rotor Rack |
| | C881-03 | Short tube adapter, for 13x75mm tubes |
| | C881-04 | Micro / Pediatric tube adapter for 0.5-2.0ml tubes |
| Specification | Power Input | 110V |
| Size and Weight | Size | 10.5 in x 12.5 in x 9 in (270 mm x 320 mm x 220 mm) |
| | Shipping Size | 14 in x 15 in x 12 in (350 mm x 375 mm x 300 mm) |
| | Weight | 15.5 lbs (C881 & C882); 16.5 lbs (C883) |
| | Shipping Weight | 19 lbs (C881 and C882) ; 20 lbs (C883) |

Additional Information