



Task Light

Gleamer[®]

Burton's all-purpose halogen light.

If your facility is searching for an all-purpose, rugged, compact halogen light, Burton's Gleamer is your solution. The Gleamer[®] was designed to meet healthcare facilities' most demanding specifications for high-intensity, cool operation, compact design and value. Offering 45,000 lux, a sturdy gooseneck, plus optional 16" (41 cm) extension arm, and superior quality components, the Gleamer provides high-intensity, maximum flexibility and extraordinary reliability. The Gleamer also features a unique venting system, energy-absorbing infrared exchanger and reflector coatings, which filter out infrared for a nearly heat-free beam in a cool, yet compact housing.




burton[®]

Features

- Easy bulb replacement
- Cool housing temperature and cool light beam; light can be in continuous use for hours without causing discomfort to physician or patient
- Convenient positioning with ergonomic positioning ridges
- High-quality, space-saving 4-caster base – safe and flexible portable movement with locking feature in 2 of 4 casters

Gleamer product dimensions

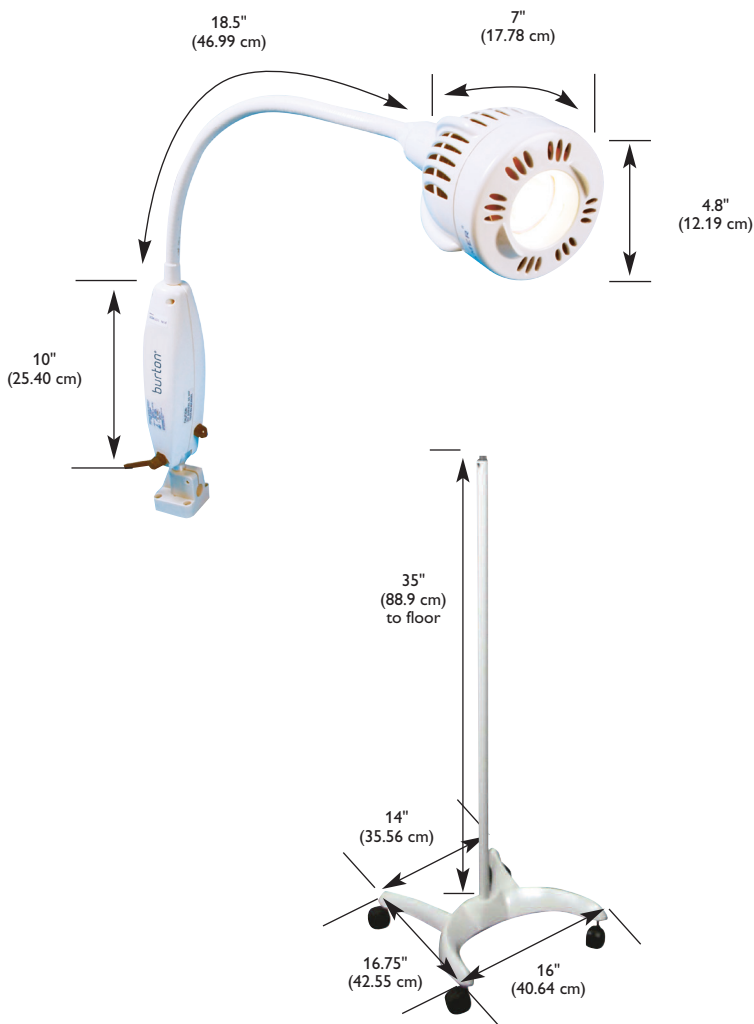

© 2018 Burton Medical. All rights reserved. Burton Medical reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Summary of technical data

GLEAMER

Illuminance:	45,000 Lux (4166 fc) at 18" (46 cm)
Color Temperature:	3700 Kelvin
Light Source:	Halogen
Rated Life:	2,000-hour average bulb life
Power:	35 Watts
Power Supply:	120v model standard; 230/240v models available
Certifications/Approvals:	UL 60601-1/IEC 60601-2-41 Certified
Product Warranty:	5-year limited warranty

Ordering Information

GLEAMER

120v Model	Description
GL30FL	Standard Spot, Floor Stand
GL30W	Standard Spot, Extension Arm & Universal Wall Mount
GLW30FL	Wide Beam Spot, Floor Stand
GLW30W	Wide Beam Spot, Extension Arm & Universal Wall Mount
230/240v Model	Description
GL32FL	Standard Spot, Floor Stand
GL32W	Standard Spot, Extension Arm & Universal Wall Mount
GLW32FL	Wide Beam Spot, Floor Stand
GLW32W	Wide Beam Spot, Extension Arm & Universal Wall Mount

ACCESSORIES

Model	Description
0007000PK	Super Bright Spot Replacement Bulbs, 4/Box

burton[®]

Burton Medical, LLC
2300 W Windsor Ct., Suite C
Addison, IL 60101 USA
Phone (800) 444-9909, Fax (800) 765-1770
Burtonmedical.com