

REACH FOR THE NEW WALLACH COLPOSTAR 1H

**A HIGH QUALITY COLPOSCOPE
THAT FOCUSES
ON AFFORDABILITY**


Take Wallach's traditional quality and dependability, add to that exceptional value, and you have the new Wallach Colpostar 1H.

The Colpostar 1H offers precision, distortion-free stereo optics, optimum single magnification, wide field of view and great depth-of-field.

Unique to single magnification colposcopes is Wallach's fiber optic halogen lighting system which yields white, clear, cool, natural tone tissue illumination for accurate colposcopic examination. In contrast, incandescent lighting used on other colposcopes tends to make tissue appear yellowish and less bright.

The Colpostar 1H features a compact center post with fine and coarse height adjustment, and a five-point base for stability and maneuverability between examination areas. In addition, Colpostar 1H meets worldwide electrical standards (UL, CSA, IEC 601).

The Wallach Colpostar 1H is the finest value in single magnification colposcopes.

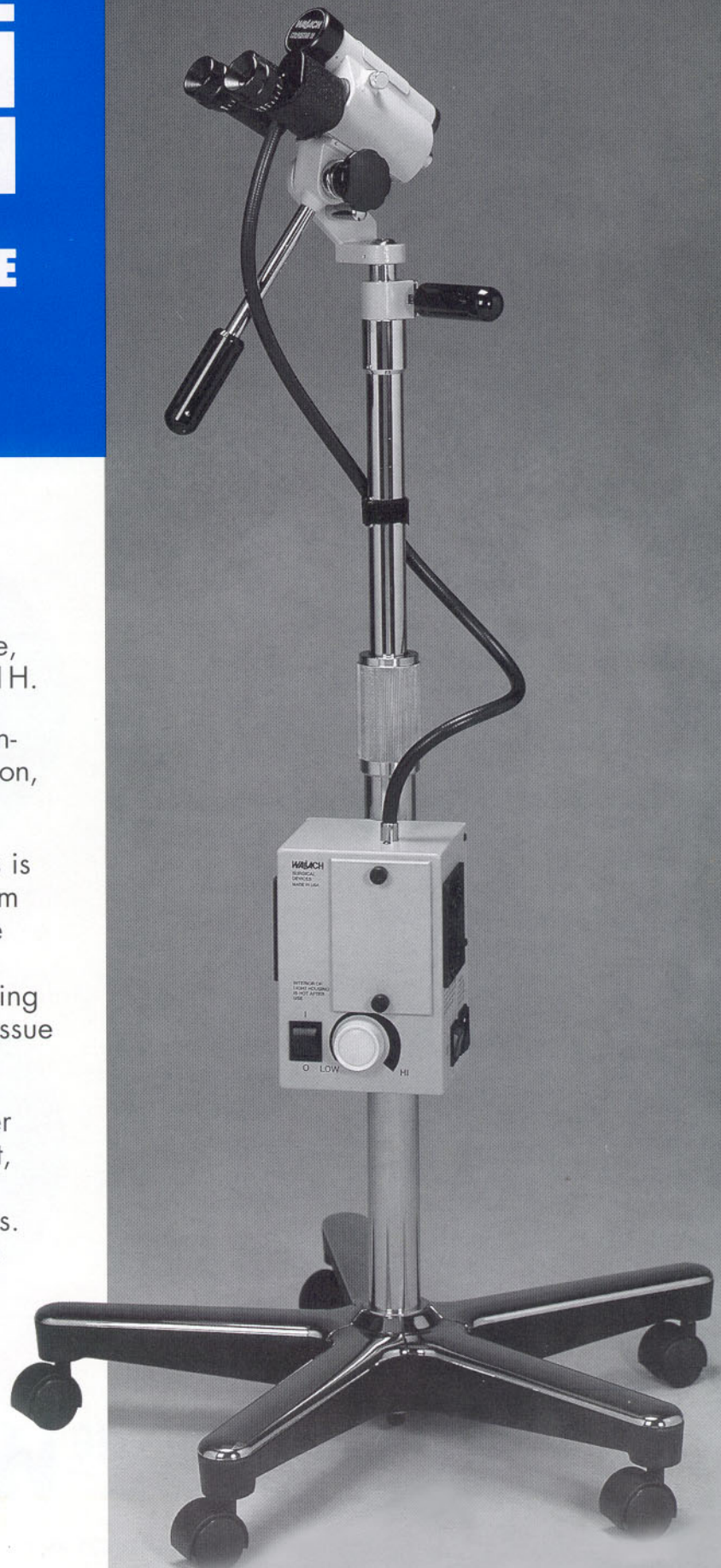

WALLACH

COLPOSTAR 1H: FEATURES & SPECIFICATIONS


1. Optical System

Binocular Eyepieces:	15x
Objective Lens:	0.8x
Magnification:	12x
Working Distance:	300mm
Visual Field:	18mm
Fine Focus Range:	40mm




2. Illumination/Power

Bulb:	150 watt Quartz Halogen
Illumination Field:	20mm
Filter:	Green Filter
Light Intensity:	1800 Foot Candles (19,368 Lux) at 300mm
Power Requirement:	2 Amps / 1 Amp 120/220 Volts AC

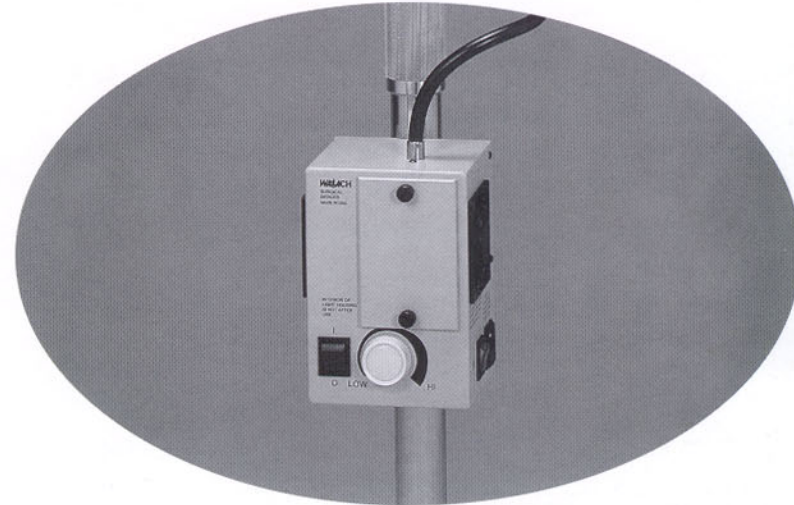

3. Stand/Base

Working Height:	86cm to 116cm
Fine Height Adjustment:	75mm
Base:	High Impact Plastic 5 Chrome covered legs with Swivel Casters

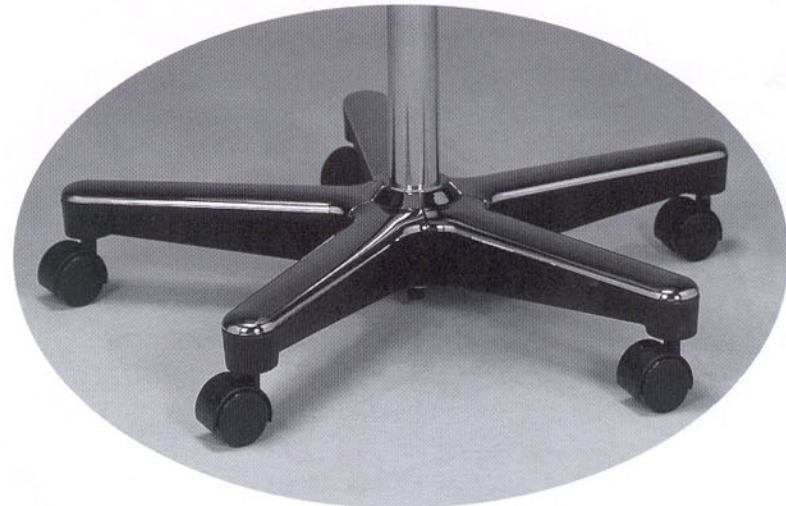

REACH FOR THE **WALLACH**