




**PRODUCT
SPECIFICATIONS**

HIGH CAPACITY PORTABLE DIGITAL SCALE

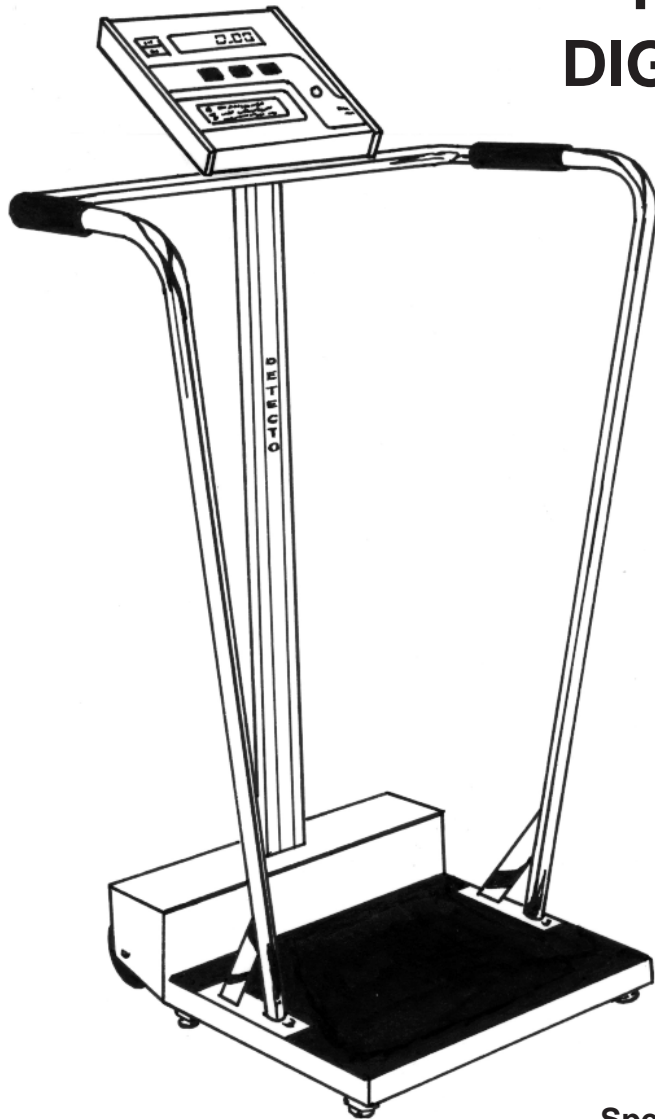

MODEL 6855

The 6855 is especially designed for obese and unsteady patients. This 600 lb capacity scale is heavy-duty, light-weight and fits through a 30" door. The handrail and large platform with its non-slip rubber mat make patient weighing quick, comfortable and secure. The LCD display can be read from either side of the platform and rotated to virtually any viewing angle. The 6855 operates on its own self-contained "C" size Alkaline batteries or with an optional power adapter for long term stationary operation.

Specifications:

Platform:	14" X 18" non-slip rubber mat
Capacity:	600 lb x .5 lb
Power:	6 "C" size Alkaline batteries (included)
Indicator:	4 digits, .7" high, 7 segment LCD
Controls:	Off/zero button
Functions:	Auto power off • Lock/release of display • Automatic zero • Motion detection
Color:	White with stainless steel handrail
Net Weight:	45 lb
Ship Weight:	70 lb
Options:	115 VAC, 60 Hz adapter with 6' power cord for long term stationary operation